


ACCENT: 800-30 & 1000-30 – Changing the Battery

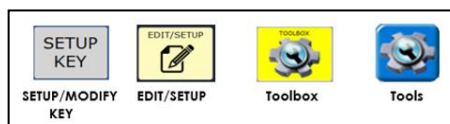

Click [HERE](#) to download NuVoice PASS

Technical Support:

Tel: +44 (0)1733 370 470 (Option 2)

Email: info@liberator.co.uk

Regional AAC Consultants:


If you wish to speak to your Regional AAC Consultant, please call +44 (0)1733 370 470 (Option 3)

The batteries in the new Accent 800-30 and Accent 1000-30 devices are now “field replaceable”, so they can now be changed by the user at home.

Removing the Device Frame

In order to change the battery, we first need to remove the frame and Keyguard/TouchGuide (if fitted) from the device.

1. If you are using a Keyguard or TouchGuide, remove it from the frame.
2. Using a Phillips screwdriver, remove the screws that secure the frame to the device (6 on an Accent 1000-30 and 4 on an Accent 800-30). **Important!** Set the screws aside. You will need them when you reattach the frame.

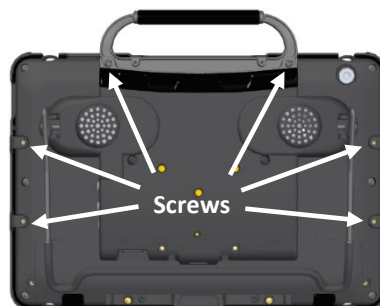

3. Grasp one side of the frame and pull it loose from the device. **Hint:** You may want to use a flat-head screwdriver or similar tool to gently lift the side and bottom tabs on the back of the device.



4. Grasp the other side of the frame and pull it loose from the device.
5. Pull the top of the frame loose.



6. Pull the bottom tab loose.



Removing the Old Battery

1. Place the device face-down on a flat surface.



2. Remove the six screws that hold the battery cover in place and remove the cover. **Important!** Set the screws aside. You will need them when you reattach the cover.

3. Grasp the plastic cable connector and pull the cable out of its socket. **Do not** pull on the wires.



4. Using the lugs to help you, slide the battery out of the device and remove from the tray. **Please contact Liberator about safe disposal of your old battery!**



Inserting the New Battery

1. Position the new battery in the tray with the cable on top and on the edge nearest you. Slide the tray back into the device.
2. Position the cable connector so that the red wires are to your right. Insert the connector in the socket. **Important!** The battery must be positioned this way. Do not try to force it to fit differently.



3. Holding the battery cover at a slight angle, insert the tabs on the cover into the slots on the back of

the device. Then lower the cover into position over the battery.



4. Secure the cover with the six screws that you previously removed.

Reattaching the Device Frame

1. Press one side of the frame onto the edge of the device.



2. Pull the other side into place. **Note:** The top tabs may snap into place as you pull on the side.



3. If necessary, press the top tabs into place.



4. Pull the bottom tab into place.



5. Look along the inside edge of the frame along the top, bottom, and sides to make sure the rubber backing didn't bunch up between the frame and the edge of the device. In this example, the rubber is bunched up, leaving a gap between the frame and the edge of the device.



If necessary, lift up the inside edge slightly and press until the rubber backing fits smoothly between the frame and device.

6. If your frame was originally secured with screws, reinsert the screws and tighten them.