


ACCENT: System Lock Settings

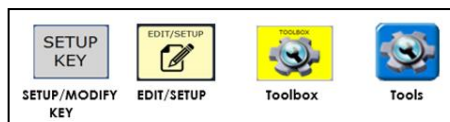

Click [HERE](#) to download NuVoice PASS

Technical Support:

Tel: +44 (0)1733 370 470 (Option 2)
Email: info@liberator.co.uk

Regional AAC Consultants:

If you wish to speak to your Regional AAC Consultant, please call +44 (0)1733 370 470 (Option 3)

When the Toolbox is locked you can still use your vocabulary to speak, but the use of the TOOLS button on the front of the case is severely limited. Only certain menus in the Toolbox will work. For example, if the system is locked and you select the TOOLS key on the front of the case, you will need to enter a password or insert the USB that has your System Lock Override tool stored on it, then select Tools.

System Lock Settings

NOTE: On older versions of the NuVoice software (v2.12 & below), and American configurations, this button is called "Disable Toolbox" and is found in "Maint. Menu".

1. Select the **TOOLS > Toolbox > Advanced Settings**.
2. Select **System Lock Settings** and turn Toolbox Restrictions **ON**:
 - a. **To create a USB lock** insert your USB stick into the device, select USB Lock, and then select Create USB Lock Override.
 - OR
 - b. **To create a Password Lock**, select Password Lock, then select Create Toolbox Password.
3. Select OK several times to exit back to your vocabulary and remove the USB stick if you have chosen this lock method.
4. The Toolbox is now locked.



If you have lost your "System Lock Override" you can download one (go to: <https://www.liberator.co.uk/support/product/system-lock-override> or select CTRL and the USB picture above). You will need to move the text file out of the zipped folder and onto the root of the USB stick (not in any other files or folders.)

