

Accent with Empower: LAMP WFL – VI Operations & Programming



 Click [HERE](#) to download Empower PASS

Technical Support:

Tel: +44 (0)1733 370 470 (Option 2)
Email: info@liberator.co.uk

Regional AAC Consultants:

If you wish to speak to your Regional AAC Consultant, please call +44 (0)1733 370 470 (Option 3)

SELECT A VOCABULARY

1. Select **Menu > Toolbox**.
 2. Select the **Vocabulary** tab, then **Select a Vocabulary**.
 3. Select the **LAMP WFL** tile. Choose your desired LAMP WFL-Visual Impairment file.
- Note:** Once you edit a button or create a Vocabulary Builder list, the vocabulary will be automatically saved and placed in **My_vocabularies**.

LAMP WFL-VI Basic Setup

You may want to change these settings to carry a black background throughout the vocabulary:

1. Select **Menu > Toolbox**.
2. Select the **Vocabulary** tab, then **Vocabulary Settings**.
3. Select **Edit Colors**.
4. Set **Home Background** and **Activity Window Background** to black.
5. Select **Button Preferences**.
6. Set Empty Button Appearance to **Empty**.

CHANGE THE VOICE

1. Select **Menu > Toolbox**.
2. From the **User Settings** tab, select **Voice**, then **Change Voice**.

EDIT A PRONUNCIATION

“Teach” the device to pronounce a word correctly.

1. Select **Menu > Toolbox**.
2. From the **User Settings** tab, select **Voice**.
3. Select **Edit a Pronunciation**.
4. Enter the correct spelling of the word, then **Finished**.
5. Enter the phonetic spelling of the word, then select **Finished**. To test the pronunciation, select the **Play** button.
6. When satisfied, select **Add**.

WORD FINDER

Use Word Finder to locate a word in the vocabulary. From the spelling keyboard:

1. Select **Word Finder**.
2. Type the word you want to find.
3. Select **Search**.
4. Select **Guide Me**. If the word is not in LAMP WFL-VI, it will say “no results found”.

SWAPPING KEYS

If an emergent communicator needs enhanced tactile cues to assist in finding highly motivating words, those words can be moved to an area corresponding to a raised location on the KeyGuide*. Once moved, those words need to remain in the new location.

1. When looking at the home screen, select **Menu > Editor**.
2. Select **Swap**. Select the buttons you want to swap.
3. Select **Save** when finished.
4. In the Full vocabulary, move the second button for the word to the corresponding location.

5. Be looking at the button you need to move. Select **Menu > Editor**.
6. Select **Swap**. Select the buttons you want to swap.
7. Select **Save** when finished.



* The KeyGuide accessory is available for Accent 1000 devices with serial numbers 5000AC10 and higher.

EDIT A BUTTON

1. Select **Menu > Editor**.
2. Navigate to the button you want to edit, then select **Edit**. Select the button you want to edit.
3. Add or edit the **Message** and the **Label**.
4. Suggested symbols will appear. If you do not see one you like, or you would like to take a **photo**, select the **Symbol** tab. You can **Search**, **Browse**, or **Import** symbols, or select **Camera**.
5. Additional button options:
 - a. **Colors:** Set various colours for the button.
 - b. **Actions:** Add an action. Action options include things such as Clear Display, Caps Lock, and Stop Speech.
 - c. **Settings:** Change button level settings. Settings include Message, Speech, Always Show Button Label, Alternate Pronunciations, and Alternate Prompt.

CREATE A BUTTON WITH A PAGE LINK

1. Select **Menu > Editor**.
2. Navigate to the button where you want a page link, then select **Edit**. Select the button you want to edit.
3. Add a **Label (in all upper case)** and a **Symbol**.
4. If desired, change the button colour from the **Colors** tab.
5. Select the **Navigation** tab.
6. If the button already has navigation, select **Remove Navigation**. Select **Go To Page**.
7. Select **Add New**.
8. **Name your page and set the grid size (84 buttons)**.
9. Select **Create**.

10. Make sure **Make Persistent** is set to **OFF**.
11. Select **Finished**.
12. Navigate to your new page and add the words you want.
See **Edit a Button**.


RESET A BUTTON TO FACTORY DEFAULTS

If you make an unintentional change or temporarily remove a page link, you can reset a button to factory defaults.

1. Select **Menu > Editor**.
2. Select **Edit** and choose the button you want to reset.
3. In the Button Edit window, select the **Settings** tab.
4. Select **Reset Button**.
5. Select **Okay** and then select **Save**.

VOCABULARY BUILDER (VB)

1. Select **Menu**, and then **Vocabulary Builder**.

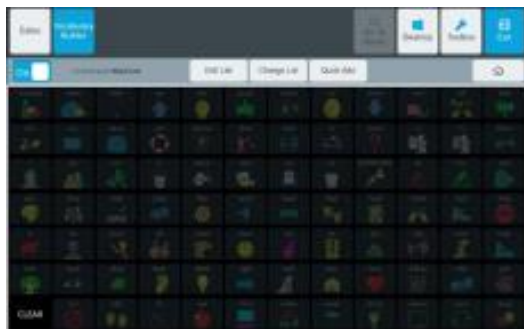


2. There are two ways to create a VB list of words.
 - a. Type a list of words.
 - i. Select **Edit List**.
 - ii. If desired, rename your list by using the text box next to **Name of Current List**.
 - iii. Select the text field and **type your word list separated by commas**.
 - iv. Select **Finished** and **Exit Menu**.
 - b. Touch the icon sequences/symbols that lead to the word.
 - i. **Touch the pathway to the first word**. Select the home button (top right) when you need to return to the home screen.
 - ii. Continue to touch the pathways to words.
 - iii. Select **Exit Menu** when finished.
3. Lists are automatically saved.

Additional VB Options:

- **Edit list:** Create, rename, manage, and clear lists.

- **Quick Add:** Use this to quickly add words to your list.
- **Change List:** If you have multiple lists, you can switch between them. In addition, rename, duplicate, merge, and delete lists.
- **My Lists and Default Lists:** In Default lists, find the LAMP WFL Transition list and activity specific lists (such as cars, bubbles etc).
- **Merge List:** If you have multiple lists, you can combine as many as you want into one list.

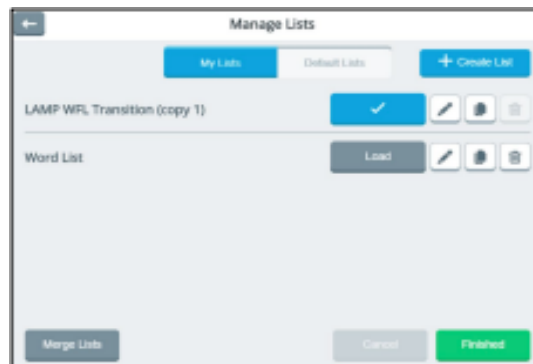

VOCABULARY BACK UP

1. Insert a USB flash drive into a USB port.
2. Select **Menu > Toolbox**.
3. Select the **System Settings** tab.
4. Select **Export & Import**.
5. Select **Export Vocabulary or Keyboard**.
6. **Export** your desired vocabulary.
7. Navigate to your USB flash drive.
8. Rename the file. It is a good idea to include the date using xx-xx-xx. For example, LAMP WFL FULL Backup 3-5-19.
9. Select **Save**.

CONNECT TO Wi-Fi®

1. Select **Menu > Toolbox**.
2. Select the **System Settings** tab.
3. Select **Devices & Connections**.
4. Select **Choose** and then connect to a network.
5. To return to Empower, select the PRC logo  that appears in the taskbar at the top.

UPDATE SOFTWARE VIA Wi-Fi

1. Select **Menu > Toolbox**.
2. Select the **System Settings** tab.
3. Select **Updates**.

4. Select **Check for Updates** next to Update Using Wi-Fi and follow the prompts.

RESTRICTIONS

To lock/unlock the Menu/Toolbox use the USB flash drive method:

1. Select **Menu > Toolbox**.
2. From the **User Settings** tab, select **Restrictions**.
3. Select **USB**. Insert a USB flash drive and follow the prompts. This process places an unlock key on the USB flash drive. Anytime this flash drive is plugged into the device, the Menu button will be unlocked.
4. **Remove the USB flash drive**, and the Menu button will be locked.

Menu Button

For LAMP WFL, keep this setting toggled **ON**.