


Accent with Empower: LAMP WFL Quick Reference Guide

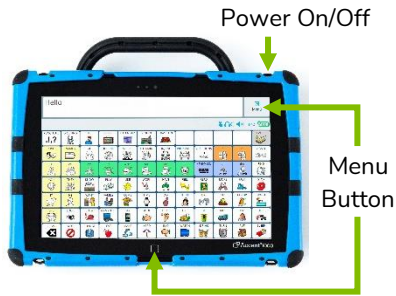

Click [HERE](#) to download Empower PASS

Technical Support:

Tel: +44 (0)1733 370 470 (Option 2)

Email: info@liberator.co.uk

Regional AAC Consultants:

If you wish to speak to your Regional AAC Consultant, please call +44 (0)1733 370 470 (Option 3)

SELECT A VOCABULARY

1. Select **Menu** > **Toolbox**.
2. Select the **Vocabulary** tab, then **Select a Vocabulary**.
3. Select the **LAMP WFL** tile. Choose your desired LAMP WFL vocabulary.

Note: Once you edit a button or create a Vocabulary Builder list, the vocabulary will be automatically saved and placed in **My Vocabulary**.

VOCABULARY SETTINGS

In this tile you can rename your vocabulary, change the font and font size on the buttons, change the background colour and set button preferences.

1. Select **Menu** > **Toolbox**.
2. Select the **Vocabulary** tab, then **Vocabulary Settings**.

RENAME A VOCABULARY

From the **Settings** tab, edit the **Vocabulary Name**.

SET BUTTON PREFERENCES (recommended settings in **bold**)

From the **Button Preferences** tab, set the following:

- **Empty Button Appearance:** What do you want empty buttons to look like? **White**, Transparent, Empty.
- **Empty Button Action:** What do you want empty buttons to do? **No Action**, Go To Home.
- **Show Labels:** Where do you want labels showing? **Home**, Middle, End.
- **Use Alternative Icons:** Where do you want an alternate set of icons to show? None, Middle, End, **All**.
- **Button Padding:** How big do you want the space between buttons to be? **0%** - 40%

CHANGE THE VOICE

1. Select **Menu** > **Toolbox**.
2. From the **User Settings** tab, select **Voice**, then **Change Voice**.

EDIT A PRONUNCIATION

“Teach” the device to pronounce a word correctly.

1. Select **Menu** > **Toolbox**.
2. From the **User Settings** tab, select **Voice**.
3. Select **Edit a Pronunciation**.
4. Enter the correct spelling of the word, then **Finished**.
5. Enter the phonetic spelling of the word, then select **Finished**.
6. To test the pronunciation, select the Play button. When satisfied, select **Add**.

WORD FINDER

Use Word Finder to locate a word in the vocabulary. From the spelling keyboard:

1. Select **Word Finder**.
2. Type the word you want to find.
3. Select **Search**.
4. Select **Guide Me**. If the word is not in LAMP WFL, it will say “no results found”.

EDIT A BUTTON

1. Select **Menu** > **Editor**.
2. Navigate to the button you want to edit, then select **Edit**. Select the button you want to edit.
3. Add or edit the **Message** and the **Label**.
4. Suggested symbols will appear. If you do not see one you like, or you would like to take a **photo**, select the **Symbol** tab. You can **Search**, **Browse**, or **Import** symbols, or select **Camera**.
5. Additional button options:
 - a. **Colors:** Set various colours for the button.
 - b. **Actions:** Add an action. Action options include things such as Clear Display, Caps Lock, and Stop Speech.
 - c. **Settings:** Change button level settings. Settings include Message, Speech, Always Show Button Label, Alternate Pronunciations, and Alternate Prompt.

CREATE A BUTTON WITH A PAGE LINK

1. Select **Menu** > **Editor**.
2. Navigate to the button where you want a page link, then select **Edit**. Select the button you want to edit.
3. Add a **Label (in all upper case)** and a **Symbol**.
4. If desired, change the button colour from the **Colors** tab.
5. Select the **Navigation** tab.
6. If the button already has navigation, select **Remove Navigation**. Select **Go To Page**.
7. Select **Add New**.
8. **Name your page and set the grid size (84 buttons)**.
9. Select **Create**.
10. Make sure **Make Persistent** is set to **OFF**.
11. Select **Finished**.
12. Navigate to your new page and add the words you want. See **Edit a Button**.

RESET A BUTTON TO FACTORY DEFAULTS


If you make an unintentional change or temporarily remove a page link, you can reset a button to factory defaults.

1. Select **Menu** > **Editor**.
2. Select Edit and choose the button you want to reset.
3. In the Button Edit window, select the Settings tab.
4. Select Reset Button.
5. Select Okay and then select Save.

VOCABULARY BUILDER (VB)

Use Vocabulary Builder to temporarily limit the number of words appearing on the screen. The strategy of “masking” vocabulary

allows a small amount of vocabulary to be taught on the device while keeping learned motor plans constant. For example, you want to start with 8 words when first introducing the device to a child.

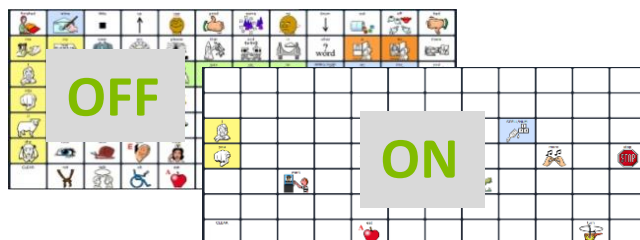


1. Select **Menu > Vocabulary Builder**. Turn VB on by selecting the **on/off** toggle.



2. Select **Hide All**. Confirm you want to hide all by selecting **Okay**.
3. There are two ways to create a VB list of words:
 - a. Type a list of words.
 - i. Select **Edit List**.
 - ii. If desired, rename your list by selecting **Change** in the upper right-hand corner.
 - iii. Select the text field and **type a list of words you want to show separated by commas**. Select **Finished** and **Okay**.
 - iv. Select **Finished** and **Exit Menu**.
 - b. Touch the icon sequences/symbols that lead to the word.
 - i. **Touch the pathway to the first word**. Select the home button (top right) when you need to return to the home screen.
 - ii. Continue to touch the pathways to words.
 - iii. Select **Exit Menu** when finished.
4. Lists are automatically saved.

Additional VB Options:

- **Edit list:** Create, rename, manage, and clear lists.
- **Quick Show:** Use this to quickly add words to your list.
- **Change List:** If you have multiple lists, you can switch between them. In addition, rename, duplicate, merge, and delete lists.
- **My Lists and Default Lists:** In Default lists, find the LAMP WFL Transition list and activity-specific lists (such as cars, bubbles etc.).
- **Merge Lists:** If you have multiple lists, you can combine as many as you want into one list.

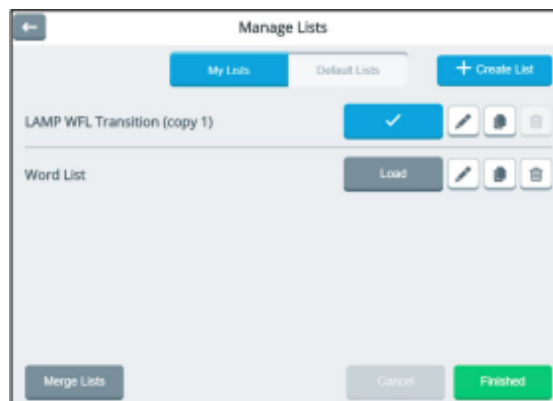

VOCABULARY BACK UP

1. Insert a USB flash drive into a USB port.
2. Select **Menu > Toolbox**.
3. Select the **System Settings** tab.
4. Select **Export & Import**.
5. Select **Export Vocabulary or Keyboard**.
6. **Export** your desired vocabulary.
7. Navigate to your USB flash drive.
8. Rename the file. It is a good idea to include the date using xx-xx-xx. For example, LAMP WFL FULL Backup 3-5-19.
9. Select **Save**.

CONNECT TO Wi-Fi®

1. Select **Menu > Toolbox**.
2. Select the **System Settings** tab.
3. Select **Devices & Connections**.
4. Select **Choose** and then connect to a network.
5. To return to Empower, select the Liberator logo  that appears in the taskbar at the top.

UPDATE SOFTWARE VIA Wi-Fi

1. Select **Menu > Toolbox**.
2. Select the **System Settings** tab.
3. Select **Updates**.
4. Select **Check for Updates** next to Update Using Wi-Fi and follow the prompts.

RESTRICTIONS

You can lock/unlock the menu two ways:

To lock/unlock via a password:

1. Select **Menu > Toolbox**.
2. From the **User Settings** tab, select **Restrictions**.
3. Select **Password**. Follow the prompts to create a password.

To lock/unlock via a USB flash drive:

1. Select **Menu > Toolbox**.
2. From the **User Settings** tab, select **Restrictions**.
3. Select **USB**. Insert a USB flash drive and follow the prompts. This process places an unlock key on the USB flash drive. Anytime this flash drive is plugged into the device, the Menu button will be unlocked.
4. **Remove the USB flash drive**, and the Menu button will be locked.

Menu Button

For LAMP WFL, keep this setting toggled **ON**.